

Elora™ Dual-Mode IPL® System for Hair Removal and Skin Treatments

Proven technology from the inventors of Intense Pulsed Light (IPL)

The Elora is a non-invasive, entry-level system for hair removal and IPL Skin Treatments using Photorejuvenation—low-risk procedures with high rates of patient satisfaction

- Simple function control, easy to learn, easy to use
- Ergonomically designed treatment head
 - Extended life: guaranteed 50,000 pulses
 - Large 20 x 30 mm spot size for rapid treatment
 - Interchangeable light guide for easy cleaning
 - Bilateral control buttons for comfortable left- or right-handed operation
- Convenient, integrated storage bin and adjustable tray
- Small and portable



Elora Dual-Mode IPL System

IPL treatments for a range of face and body skin conditions

- Vascular changes
 - Redness
- Pigmentary changes
 - Sun damage
 - Age spots (sun-induced freckles)
 - Brown spots
- Hair removal for all body areas and hair depths



Elora™ Dual-Mode IPL® System

Specifications

Light Source	IPL (Intense Pulsed Light)
Wavelength Range	570-1200 nm
Applications	Hair removal / IPL Skin Treatments using Photorejuvenation
Treatment Parameters	3 predefined programs
Energy Density (fluence)	Up to 22 J/cm ²
Spot Size	20 x 30 mm
Pulse Repetition Rate	0.33 Hz
Dimensions	38 x 38 x 60 cm / 15 x 15 x 24 in (WxDxH)
Weight	20 kg / 44 lbs
Electrical Requirements	95-252 VAC, 5A Single phase, 50/60 Hz


The Elora intense pulsed light system treats a broad range of applications safely and effectively. This small, portable and simple-to-operate system readily integrates into any practice.




Before lentiginos on the hand



After four IPL treatments using the Elora

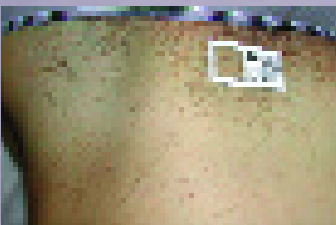

Before bikini line and thigh hair

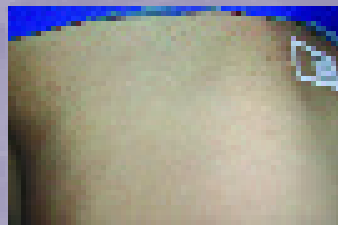

After five IPL treatments using the Elora

Lumenis Global Sales Offices

Lumenis Inc.
2400 Condensa Street
Santa Clara, CA 95051 USA
Tel 408 764 3000
877 586 3647
Fax 408 764 3999
800 505 1133

Lumenis (GERMANY) GmbH
Dreieich-Dreieichenhain,
Germany
Tel + 49 6103 8335 0
Fax + 49 6103 8335 300

Lumenis (FRANCE) SARL
Orsay, France
Tel + 33 1 69 33 14 20
Fax + 33 1 60 19 57 23

Lumenis (ITALY) Srl.
Formello (RM), Italy
Tel + 39 06 90 75 230
Fax + 39 06 90 75 269

Lumenis (UK) Ltd.
London, UK
Tel + 44 20 8324 4200
Fax + 44 20 8324 4222

Lumenis (JAPAN) Co. Ltd.
Tokyo, Japan
Tel + 81 3 5789 8300
Fax + 81 3 5789 8310

Lumenis (CHINA) Ltd.
Beijing, China
Tel + 86 10 6510 2620
Fax + 86 10 6510 2621

Lumenis (HK) Ltd.
Hong Kong
Tel + 852 2174 2800
Fax + 852 2722 5151

Internet: www.lumenis.com

Email: information@lumenis.com