


□ one platform □ □ multiple world-class technologies



SapphireCool™ Light Guides for Universal IPL® Treatment Head

15 x 35 mm

8 x 15 mm

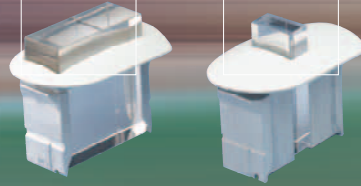

Single Universal IPL® Treatment Head for all IPL Applications



lumenis **1one**™
Universal IPL®

Technological Advances with the Universal IPL® Treatment Head

- **From the inventors of patented IPL® technology**
 - 5-year, long-term clinical studies available on safety and efficacy of Lumenis IPL
- **Revolutionary Universal IPL® treatment head**
 - One single, versatile, IPL head for all IPL applications
 - Easy to use; no need to change treatment heads when changing applications
 - Space efficient design eliminates need to store multiple treatment heads
- **7 ExpertFilters™ enable optimized treatments**
 - 515, 560, 590, 615, 640, 695, 755 nm filters
 - Computer-aided filter recognition
- **Built-in calibration**
- **Up to 3.5 times faster than other systems**
 - 1 second repeat, large spot size (15 x 35 mm)
 - Results in greater patient throughput
- **Continuous contact cooling for increased patient comfort**
- **Low cost of ownership**
- **2 SapphireCool™ Light Guides**
 - Increased patient comfort
 - Enhanced durability of the light guide
 - 15 x 35 mm and 8 x 15 mm spot sizes for efficient treatment at every treatment site
- **Ergonomic design reduces user fatigue**
 - 3 operating buttons for multiple gripping positions
 - Equally efficient for left and right handed users
- **Dossi™ Patient Database**
 - Efficient storage of patient information and treatment parameters
- **A multitude of applications**
 - IPL® skin treatments using Photorejuvenation, vascular and pigmented lesions and hair removal
- **Easy to learn, easy to use**
 - Pre-set parameters in Basic and Advanced modes
 - Advanced mode for complete customization of treatment parameters

lumenis **one**
Universal IPL®

Lumenis Global Sales Offices

Lumenis Inc.
2400 Condensa Street
Santa Clara, CA 95051 USA
Tel 408 764 3000
877 586 3647
Fax 408 764 3999
800 505 1133

Lumenis (GERMANY) GmbH
Dreieich-Dreieichenhain, Germany
Tel + 49 6103 8335 0
Fax + 49 6103 8335 300

Lumenis (FRANCE) SARL
Orsay, France
Tel + 33 1 69 33 14 20
Fax + 33 1 60 19 57 23

Lumenis (ITALY) Srl.
Formello (RM), Italy
Tel + 39 06 90 75 230
Fax + 39 06 90 75 269

Lumenis (UK) Ltd.
London, UK
Tel + 44 20 8324 4200
Fax + 44 20 8324 4222

Lumenis (JAPAN) Co. Ltd.
Tokyo, Japan
Tel + 81 3 5789 8300
Fax + 81 3 5789 8310

Lumenis (CHINA) Ltd.
Beijing, China
Tel + 86 10 6510 2620
Fax + 86 10 6510 2621

Lumenis (HK) Ltd.
Hong Kong
Tel + 852 2174 2800
Fax + 852 2722 5151

internet: www.lumenis.com
email: information@lumenis.com



Revolutionary Universal IPL® Treatment Head



Advanced Mode Treatment Screen for IPL® Skin Treatments

Clinical Evidence

IPL Skin Treatments on Face



Before
Basmat Zamir, MD, Haifa, Israel



After 2 treatments

IPL Skin Treatments on Upper Chest



Before
Basmat Zamir, MD, Haifa, Israel



After 1 treatment

Technical Specifications

Applications	IPL Skin Treatments; Vascular Lesions; Pigmented Lesions; Hair Removal
Spectrum	515-1200 nm
Energy Density (Fluence)	10-40 J/cm ²
Pulse Duration	3-100 ms
Pulse Delay	1-120 ms
Pulse Characteristic	Multiple Synchronized Pulsing
Repetition Rate	Up to 1 Hz
Spot Size	15 x 35 mm, 8 x 15 mm
Cooling	Continuous contact cooling

System Specifications

Size	27" x 18" x 63" (67 x 47 x 159 cm)
Weight	338 lbs (154 kg)
Electrical Requirements	100-230 VAC 10%, 25/16A, 50/60 Hz Single Phase. Dedicated Line.



www.aesthetic.lumenis.com

Caution: Federal law restricts this device to sale by or on the order of a physician.

Lumenis and its logo, Lumenis One and its logo, Aluma, LightSheer, IPL, MultiSpot, Dossi, SapphireCool, ExpertFilter and ChillTip are trademarks of the Lumenis Group of Companies. Product specifications are subject to change without notice.

Copyright © 2005, the Lumenis Group of Companies. All rights reserved. PBP000321 Rev C



0473